



CONTIPHARMA


a member of Bedelco group

PROTECT

Disposable surgical mask Purism Type IIR



OUTSIDE

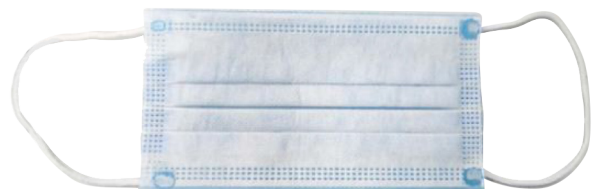

INSIDE

BRAND

PURISM

BACTERIAL FILTRATION EFFICIENCY

≥ 99%



SPECIFICATIONS

3-layer mask :

- First layer: waterproof non-woven fabric (new non-woven fabric): effectively inhibits and insulates harmful particles and gases
- Second layer: high-density filtering layer (blown non-woven fabric)
- Third layer: layer in direct contact with the skin, in new non-woven fabric, resistant to moisture, allowing easy, non-toxic, non-stimulating, soft and skin-friendly breathing

Composition :

unwoven-bound fabric, fusion-blown non-woven fabric, earrings and nasal bridge bar.

Recommended port time :

4-8 hours, to be changed between medical procedures. After contact with the used mask, be sure to disinfect your hands.

APPLICATION

for medical personnel during operations. It covers the user's mouth, nose and jaw, providing a physical barrier to prevent the direct passage of pathogens, microorganisms, bodily fluids and particles. Don't use in a traumatic environment.

COMPLIANCE

EC marking according to European Directive 89/686/EEC
Standard EN 14683 Type IIR

CONDITIONING

10 pcs/ bag, 5 bags / box, 48 boxes / cardboard, 2400 pieces / cardboard,
Cardboard size: 47.5 × 41.5 × 56 cm - 11.7 kg

STORAGE

to be stored in a dry, clean, ventilated and non-corrosive gaseous area, away from ignition sources or combustible materials

STORAGE TIME

3 years

